PRIMAX[®] RTG-B

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Density curve:
see diagram


Average gradient:
3,2 (+/- 15%) as per DIN 6867, part 1

Darkroom light:
Red safety filter with light intensity of 57 mlux (DIN 6868/part 2)
Filter recommendation:
Agfa R-1, Kodak GBX-2 or similar ones

Processing:

- Standard 90 seconds processing system at 34°C with minimum 90 seconds developing/fixing/washing/drying or according to recommendation of the producer
- Manual processing at 20°C/5 minutes developing time (hardener required to be added to the fixing bath). If using other temperature/time combination process must be adopted.

In case of manual processing it is to be considered that contrast, speed and maximum density are differing from the figures given herewith for the machine processing.

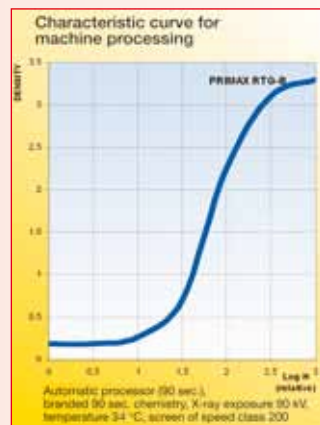Screen recommendation:

- Calcium tungstate:
PRIMAX B 200,
CAWO Universal, Special
Agfa CPB 100, CPB 200
- Rare earth:
PRIMAX Y 400
CAWO SE 4, SE 8
Agfa Curix Blue 200, 400, 800
Agfa CPB 400
Kodak X-omatic Regular

Packing and Sizes:

- All usual cm and inch sizes are available
- Packing in 100 sheet boxes, 5 boxes of 100 sheets in one transport master carton
- Packing NIF (without interleaves between film sheets)

Storage:

- unprocessed film should be stored at maximum 23°C
- relative humidity to be kept within 40 to 60 per cent
- protection necessary against ionizing radiation (max. 90 nGy/h)
- protection necessary against the influence of chemicals

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